

DETECT™

APEX LOCATOR

NEW GENERATION APEX LOCATOR WITH FULL COLOR DISPLAY



SHARPEN YOUR ROOT CANAL VISIBILITY

DENTSPLY
MAILLEFER

DETECT™

APEX LOCATOR




PRODUCTS FEATURES

- Large color display to facilitate file tracking
- Progressive sound control for dual control
- Powered by rechargeable battery
- Latest multi-frequency technology enables location of the apex in most types of canal conditions
- Fully automatic device


UPPER CANAL ZONE

The cursor on the tooth icon indicates the progression of the file in the canal

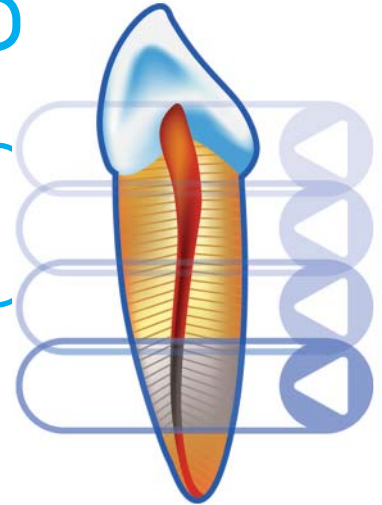


APICAL ZONE


Extended apical zone from 0.9 to 0.0

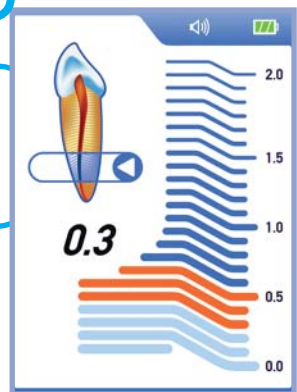

When the apex is reached, a solid tone is emitted

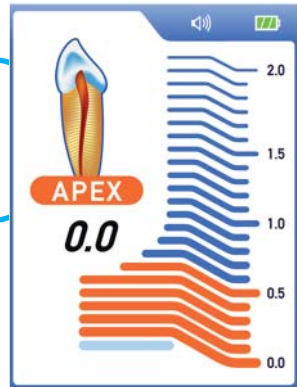

When the cursor on the tooth icon shows **OVER** this indicates that the file has passed the apex

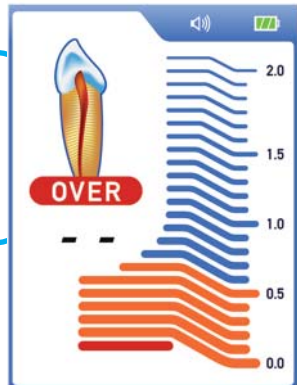

+
WE
KNOW
ENDO.